

# Foot Massager Machine

## Instruction Manual

Thank you for purchasing our products.

- Please read the instruction manual carefully before use, especially the "safety precautions", fully understand the scope of application and operation methods of the product, so that you can use it correctly and safely.



- Persons who have recently undergone surgery
- Persons with sensitive skin or poor blood circulation

If you have any of the following descriptions, you must consult your doctor and ask if you can use it:


- Persons with malignant tumors (tumors)
- Persons with heart disorders
- For those who have lost temperature sense
- Persons with sensory impairments
- Persons with wounds at the treatment site
- People with a body temperature of 38°C or higher (e.g., acute inflammation symptoms, fatigue, chills, blood pressure fluctuations) people who need to rest or when their physical condition is significantly poor.

## Product Parameter




Product Model	HZ-505
Product Name	Foot Massager Machine
Timing Function	5min \ 10min \ 15min \ 20min \ 25min \ 30min
Heating Temperature	40-45°C
Vibration Function	F1
Comprehensive Mode	P1 \ P2 \ P3 \ P4 \ P5
Air Pressure Intensity	L1 \ L2 \ L3 \ L4 \ L5 \ L6 \ L7 \ L8 \ L9
Function Mode	Foot Sole Kneading \ Air Compression Massage \ Plantar Vibration \ Plantar heating
Rated Power	50W
Rated Voltage	110V (US) / 220-240V (EU/UK)




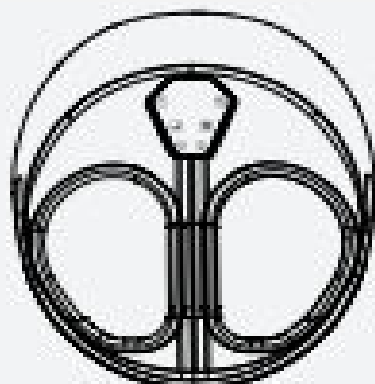

## Product Description

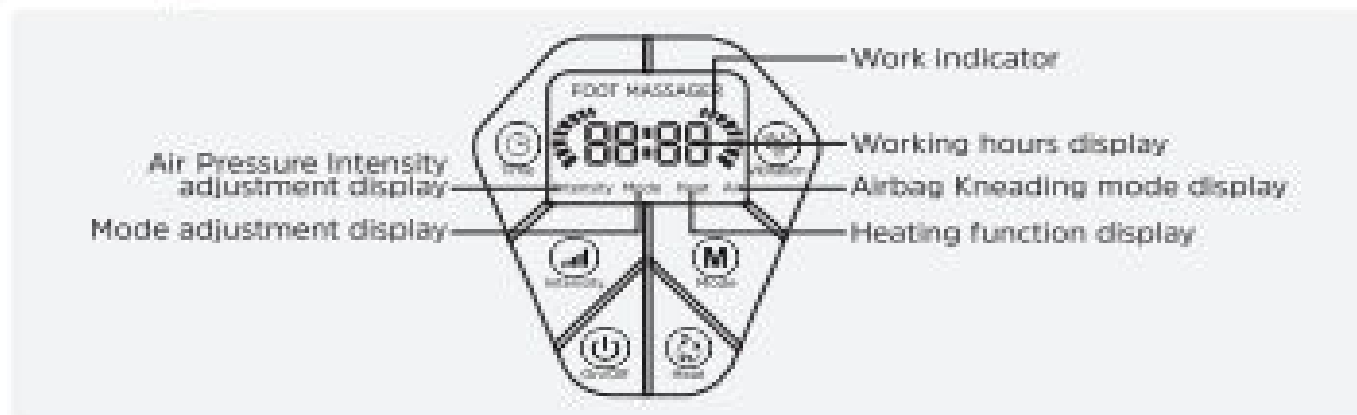

## Screen Display Introduction:

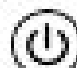

## Introduction of button function:

-  : Control product start-up and shutdown.
-  : The switch to control the heating function, short press to turn on the heating function, the screen displays "heat"; short press to turn off the heating function, no display on the screen.
-  : Control the massage time of each massage mode. Short press to select the time. There are 6 time periods of 5min, 10min, 15min, 20min, 25min, and 30min to choose from. The screen will display the selected massage time. After the selection, it starts to enter the use time countdown.
-  : Press to turn on the vibration function, the vibration mode is not turned on in the "FO" state, short press to turn on the vibration mode, and at the same time, "FO" on the screen will change to "F1".
-  : Choose the massage mode, there are 5 modes for you to choose.

P1	Airbag Kneading working + roller rotates and massages the feet in 60 seconds with a slow frequency and intermittently holding the feet, roller switch a forward and reverse cycle every 60 seconds.
P2	Airbag Kneading working + roller rotates and massages the feet in 30 seconds with a middle frequency and intermittently holding the feet, roller switch a forward and reverse cycle every 30 seconds.
P3	Airbag Kneading working + roller rotates and massages the feet in 10 seconds with a middle frequency and intermittently holding the feet, roller switch a forward and reverse cycle every 10 seconds.
P4	Airbag Kneading off + roller rotates and massages the feet in 60 seconds with a slow frequency, roller switch a forward and reverse cycle every 60 seconds.
P5	Airbag Kneading working + roller stop working.

### Note:

In P4 mode, only the rollers run, the airbag is turned off, and "air" is not displayed on the screen.

-  : Select the pressure of the airbag. There are 9 pressure levels to choose from L1 to L9, the pressure gradually increases.

Massage intensity (9 levels)	L1	L2	L3	L4	L5	L6	L7	L8	L9
------------------------------	----	----	----	----	----	----	----	----	----

**Note:**

- The default setting is L5 when the product is turned on, and the massage intensity can be adjusted according to your needs;
- L1-L3 is suitable for the elderly who like gentle massage; L3-L6 is moderate in strength, suitable for most people; L7-L9 is stronger in strength, suitable for people who need to massage acupoints with more force.

**Introduction of wireless remote control:**











- When you use it, it is inconvenient to bend over, you can directly use the wireless remote control to operate.
- The button function of the wireless remote control corresponds to the button function of the main unit. The wireless remote control does not have a separate button to control the vibration function. Please operate on the main product to turn on and off the vibration function.

**Note:**

- The wireless remote control should be used within 1.5m.
- You need to purchase and install the AAA battery by yourself.

## Operation Guide



1. Please check whether the product is damaged or not before use, and then insert the plug into the power socket after confirming that it is placed in a stable manner.
2. When you see a blue light bar on the display of the machine  , it means that the device is already powered on.
3. After digging into the feet to adjust the appropriate position, short press  to turn on the device; at this time, the default massage mode is P1, the massage time is 15 minutes, and the airbag strength is L5.
4. During use, you can press  to select the desired massage mode; press  to select the appropriate airbag compression strength; press  to turn on or off the vibration mode; press  to turn on or off the sole heating mode; through press  to select the massage time.
5. When the massage time is up, the product will automatically shut down. If you want to shut down in the middle, you can directly press  to shut down the device.

