


MASSAGE GUN

USER MANUAL



FUNCTION OF KEYS

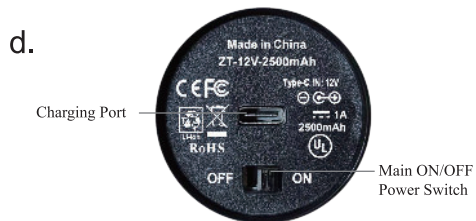

INSTALLATION OF MASSAGE HEADS



Pull it out do not rotate.



Insert the massage head completely by pressing.


MASSAGE HEADS



- **Round Head**
Relax large muscle group
Back, leg and waist



- **Spade Head**
Relax Hamstring, glutes,
Upgrade flexibility

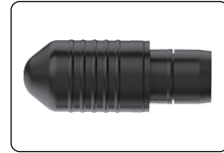

- **Long Taper**
Relieve deep soft tissue
Pain, knees and joints



- **Flat Head**
Relax and reline whole
Body muscles



- **U-Shape Head**
Relieve cervical spine
stiffness



- **Semicircle Head**
Relieve Should and arm
soreness



• **T shaped Head**

Used for arms,foot massage



• **Mushroom Head**

Suitable for relaxing joints
and sensitive muscles



• **Gyro Head**

Suitable for massging
soft tissue



• **Fan-shaped Head**

Suitable for back muscles



• **Dispersed shaped Head**

Suitable for shoulders massage



• **Thumb Head**

Uesd for massage soft tissue

OPERATING INSTRUCTIONS

1. Choose a suitable massage head based on your need (Page 3) and insert into the round hole in front of the massager (Page 2)
2. Turn on the main power switch located at the bottom of the handle (Page 1d).
3. Single tap the On/Off button (Page 1b) located on the LCD screen.
4. Continue to click the switch button to adjust to the appropriate gear.
5. After massage, please remember to turn off the main switch at the bottom of the handle.

CHARGING INSTRUCTIONS

1. Before using it for the first time, please charge it fully. The charging time usually lasts about 3 to 4 hours.
2. During the charging process, do not use the device to avoid potential damage to the battery.
3. Please charge the massage gun before its battery level drops below 20% to extend the battery life.
4. If not used for a long time, the built-in battery of the product should be charged regularly (About once a month).
5. Do not charge the device in a high temperature humidity environment or under direct sunlight.

HIGHLIGHT

• Timer Function:

If you use the massage gun continuously for 10 minutes, the device will automatically pause to avoid overheating of the motor and battery. If you want to continue using it, after the device has cooled down, you can turn on the switch again and continue using it.

• 12 Massage Heads & 7 Speed Levels:

The massage gun (model: Y12) is equipped with 12 massage heads allow you to massage every muscle tissue, which will provide you a different relaxing massage experience. Also it comes with 7 speed levels, you can fine tune the speed to achieve the most comfort and effect at different spot with different level of pain. It's not only suitable for warm-up at low speed level, but also suitable for deep tissue relaxation at high speed level.

SPECITICATION

| | |
|-----------------------------------|----------------------|
| Model: Y12 | |
| Net Weight: 0.9 kg | Gross Weight: 1.71kg |
| Adapter: 12V/1.0 A | Battery: 2500 mAh |
| Control: LCD Touch Screen | Massage Heads: 12 |
| Work Time: up to 30 cycles (10 ') | Charge Time: 2-3 H |

SAFETY TIPS

1. Children please use under adult supervision, pregnant women and the elderly use with caution and follow doctor's advice.
2. Choose the right massage head based on the body parts as described in page 3 and page 4.
3. This device should be gently pressed and moved. The massage time of the same body area should not exceed for 1 minute.
4. High vibrating speed should only be applied to soft tissue parts of the body, use with caution. No high speed and high pressure on head and bones.
5. If you feel painful and uncomfortable with the massager, stop using it immediately.
6. Avoid direct contact with the telescopic rod and heat dissipation hole of the massage gun to prevent any potential damage.

WARRANTY

Our product has 6-month warranty from the date of purchase if it fails due to a defect in mechanics or material. We provide free repair and replacement based on the damage.

Notably, no warranty for the following situations:

1. Damage to this equipment caused by incorrect operation or storage;
2. Unauthorized removal or repair of this equipment;
3. Unable to follow the instruction manual;
4. Abuse of product;

Thank you for purchasing our product. Hope you enjoy it!

Important:

1. Please do not use the product while it is charging.
2. Before charging, please turn off the main power switch.
3. Do not use or charge the device in a high temperature environment.
4. Please charge the device in time before the battery power of the device drops below 20%.

Remark:

1. Please do not dismantle the device.
2. Please use the original charger and the attachments.
3. If the device rattles or the massage head comes loose, please switch it off immediately.



Maintenance & Storage:

1. Use slightly dampened towel to clean the device surface and a soft tissue to make it dry.
2. Each time after massage, remove the massage head and put them back to the case.

COMMON PROBLEM

1. Why can't the fascia gun start?

- ① There is an on / off switch at the bottom of the battery to ensure that the on key is on.
- ② It may be unable to start when the battery is low. Try again after charging.

If the fascia gun still cannot be started, please contact customer service or sales personnel

2. What is the reason why the sound becomes louder after using it for a period of time?

The louder sound is related to the use time and frequency. If it is often used for a long time, the sound may become louder after a period of time, which is a normal phenomenon. If there is obvious abnormal noise at the motor, and you feel that the screw falls off or the support is broken, please contact the customer service or sales personnel.

3. What is the function of the light above the rear cover?

Above the back cover is a pressure sensing display lamp, which is used to reflect the strength. Generally, if no pressure is applied, the lamp above the back cover does not light up. When pressure is applied, the pressure sensing display lamp above starts to light up as the strength increases. The greater the strength, the more lights will light up, A total of 7 pressure sensing display lights help to remind you not to exert excessive pressure on your body and better grasp the strength to avoid damage.

4 . What is intelligent enhancement?

The function of intelligent strength enhancement is that in the process of use, the greater the resistance of the motor, the motor intelligent control board Instantly increases the number of revolutions to enhance the strength, so that users can feel more powerful power.

