


Shiatsu Multi-Function Massage Cushion



Shiatsu Massage Cushion

+ Structure

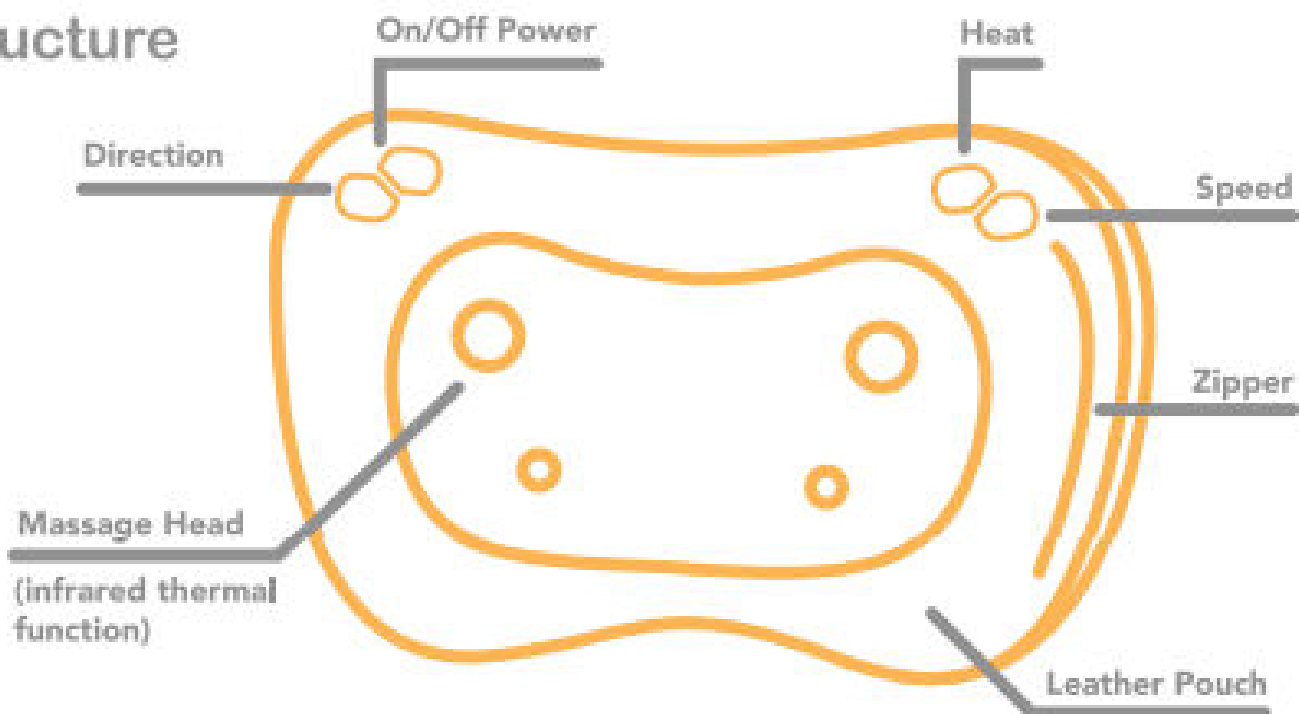

includes AC Adapter & Car power adapter*

Shiatsu Massage Cushion

+ Application Environment

- In order to avoid electric shock or malfunction, do not use the product in damp or dusty environments.
- Do not use the product under the temperature high than 40°C
- Do not use the product outdoors

+ Non-Applicable Population

If you are dealing with the following situations, consult a doctor before use.

- You are receiving medical treatment or you feel uncomfortable.
- Cancer patients
- Those who are suffering from acute diseases.
- Those who have heart pains.
- Osteoporosis patients.
- Those who are using cardiac pacemakers and or other implanted electronic medical devices

A2

Shiatsu Massage Cushion

+ Product's features & effiacy

The full body Shiatsu massage cushion, is an in-depth massage from neck, back, shoulders, legs, and waist. The heating feature can relax your body, eliminate fatigue, relieve pressure and promote blood circulation.

Compact, stylish design, not only to meet your household demand but also take into account the needs of the driver, so that you can enjoy the relaxation of the massage in the car to achieve the perfect massage effect for the body.

Shiatsu Massage Cushion

+ Notes before usage

- Do not use with other therapy equipment or electric blanket.
- Do not use within half an hour after the meal.
- To avoid injury, do not allow children use.
- Avoid carrying heavy objects with the product, causing mechanical failure or damage due to fallen objects.
- This product is designed for single person use. Do not simultaneously use with two or more people.
- This product is designed for people who weight less than 120KG.

+ Cases for ending the usage of this product

In the following circumstances you should immediately stop using this product, turn off and unplug from socket.

- Accidentally splashed water on massage cushion and or adapter.
- Related accessories are damaged or internal parts exposed.
- You feel unusual and or severe pain during or after usage.
- Any abnormal situation occurs when the massage is working.
- A sudden power failure.
- Lightning occurs.

A3

Shiatsu Massage Cushion

+ Product's features & effiacy

Description	Shiatsu Massage Cushion
Normal Operating Voltage	DC 12 V
Power Consumption	24 W
Supply Unit Data	
Description	Supply Adapter
Normal Operating Voltage	110~230 v/ 50~60 Hz
Normal Output Voltage	DC 12 V/ 2000 mA
Shock Protection Class	Class II

We hope you enjoy your new massager.

For assistance, email: service@alphamarts.com, or call our after sales support 877-606-2818,
Monday-Sunday, 9am - 6pm PST

Shiatsu Massage Cushion

+ Assembly & Repair

- Do not dismantle the sheath of the product and contact internal components of the product to avoid personal injury due to electric shock or machine failure.
- This product is only repaired by professional service personnel. Do not dismantle or repair it by yourself.

+ Notes for power plug & power wire

- Do not pull off adapter with wet hands to avoid electric shock hazard.
- Check whether it meets specified voltage for the product.
- When you pull off adapter, you should grasp the whole adapter instead of the wire.
- Timing function can automatically stop the massage cushion, please close the power switch after the use.
- If adapter is damaged, immediately stop using the product, turn off, pull out adapter from socket, and contact qualified service personnel.
- When not in use for extended period of time, turn off the power switch and power wire.

Shiatsu Massage Cushion

+ How to use

1. Insert the home adapter to connecting head and insert the power plug to the socket with 220V/50hz source, it is the power-on standby when "OFF" is indicated on the digital display
2. When use in the car, connect to the car cigarette lighter dedicated connection line, insert to the 12V power source, it is in power-on standby position when "OFF" is indicated on the digital display.
3. Press power switch key, massage function on, pressed again for heat to shut off. Cycle selection. (Pay attention to temp when using this function.)
4. Press power switch key for 3 seconds, all functions will shut off.
5. The product can switch rotation direction (clockwise rotation) automatically.
6. This product will enter the 20 minute countdown state after activated.
7. After use, please unplug the product and store it.

+ Notes

- 1. After usage, cut off the power supply.
- 2. Do not clean it with water or other corrosive chemical agents. During daily maintenance, wipe clean with soft cotton fabric.
- 3. Should avoid contact with sharp and or hard objects. Do not press or step on.
- 4. When using, if there is any abnormal situation, cut off the power. It is forbidden to modify the power adapter or repair it at discretion after the use.