




User's Manual

Neck massager

Product Components

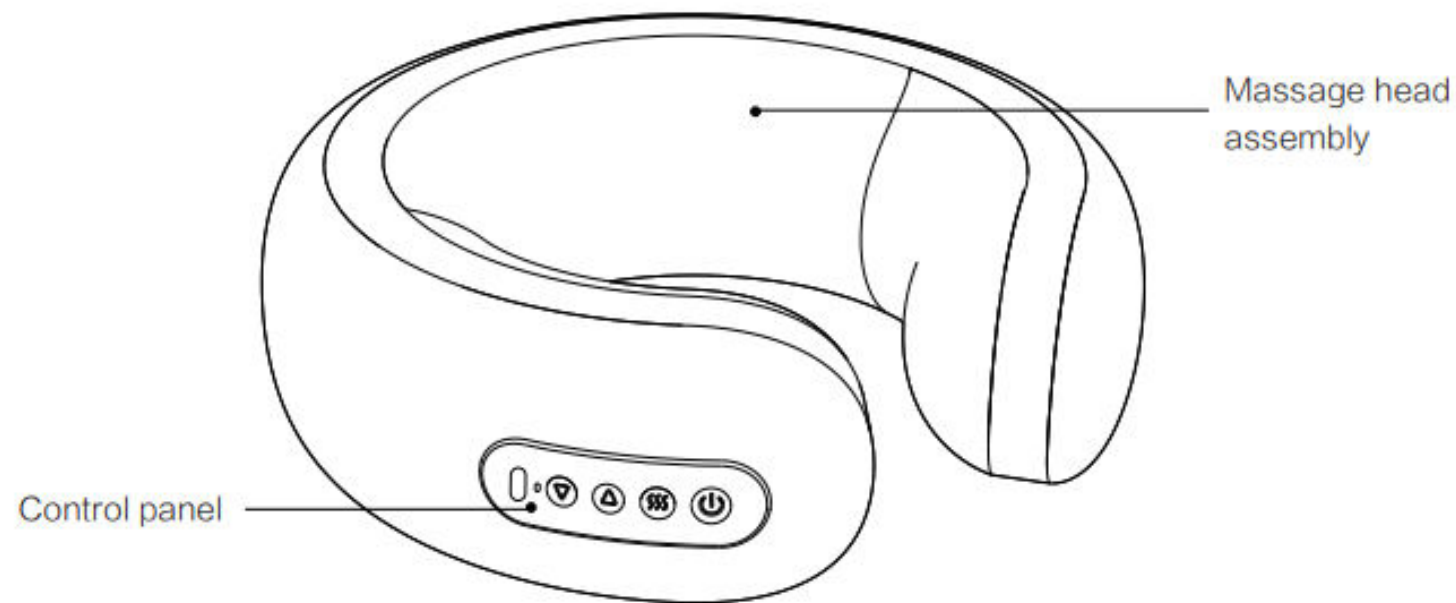

Product Information

Product Name: U-Shape Massage Pillow
Rated time: 15min
Rated voltage: 7.4V ---
Host size: 27x25.8x9.5cm

Product model: HT-PN031
Battery capacity: 2200mAH
Rated power: 6W
Package size: 26x25x9cm

Applicable Standard: GB4706.1-2005, GB4706.10-2008


The product is subject to the actual product. Thank you for purchasing this product. Please read this manual carefully before use.

Charging

When the massager is short of power, one end of Type-C cable can be inserted into the massager, and the other end can be inserted into the phone charger with Android head (the output is 5V --- 2A), red on charging, green after full charge.

- The red light would stay on throughout the charging process.
- The green light would stay on if the product is fully charged.
- The green light would stay on if the power is more than 20%.
- The red light would flash if the power is less than 20%.

Function of the Control Pane

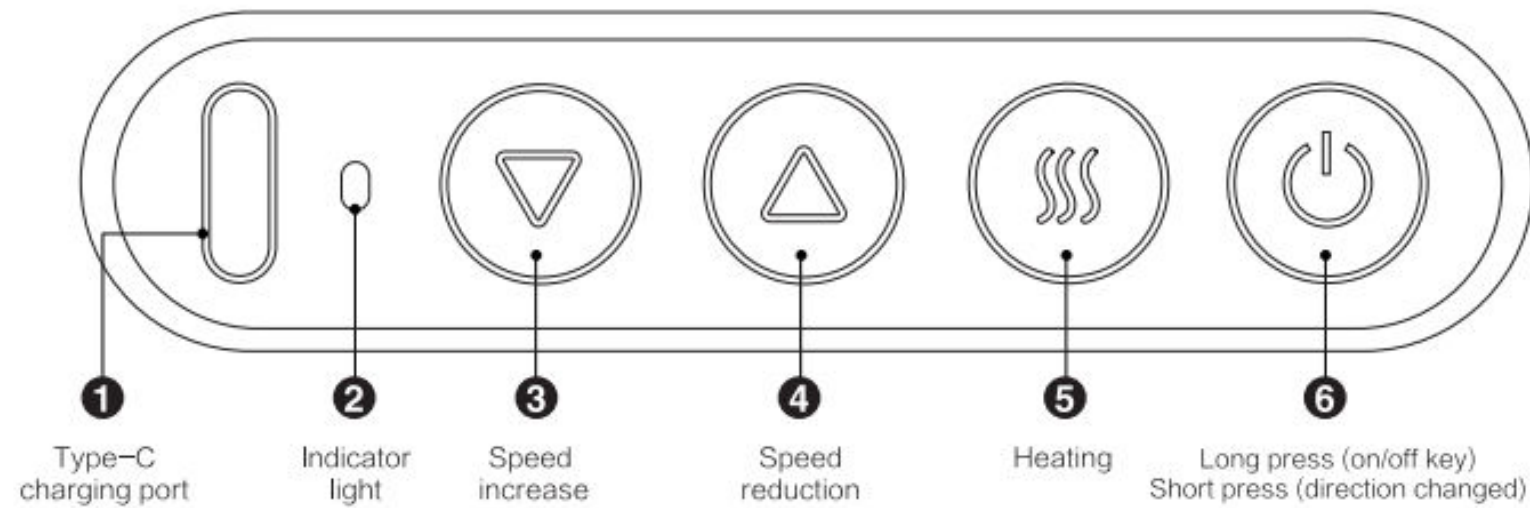

Operation method



1. Power key

Long press the power button to (start/stop) the massage function. Short press the button to adjust the forward and reverse massage.



2. Heating key

Short press the heating button after starting up, (start/stop) the heating function.



3. Gear key (speed reduction)

Short press the gear button after turning on the machine to reduce the massage speed.



2. Gear key (speed increase)

Short press the gear button after turning on the machine to increase the massage speed.

Attention:

- > It is recommended that the massage time not exceed 15 minutes.
- > Please power off the massager if you feel uncomfortable during the massaging process.
- > Please power off the massager if you do not need it.
- > If any malfunction occurs, please turn off the main switch to avoid accidents.

Precautions

! Dangers

Please never use it with the following medical electronic devices:

1. Implantable medical electronic devices such as cardiac pacemakers;
2. Life-sustaining medical electronic devices such as artificial hearts and lungs;
3. Electrocardiograph and other medical electronic equipment;
If used together, these medical electronic devices will cause malfunction, which may be life-threatening;
4. Gas stations, high-dust workshops, gas and other flammable and explosive places are prohibited;

! Warnings

If you feel unwell during or after use, stop using it and consult a doctor. This product contains car batteries, please do not disassemble, modify, immerse the product in water or throw it into fire, otherwise the product may be damaged and explode.

Please consult your doctor if any of the following circumstances, otherwise, incidents may happen or you may suffer from physical discomfort:

1. Patients who are receiving medical treatment and feel abnormal in their bodies.
2. Patients with heart and cranial nerve abnormalities.
3. Patients whose body temperature exceeds 38 degrees (fever period).
4. Patients with infectious diseases.
5. Patients with impaired skin perception or skin abnormalities.
6. Patients who are taking drugs.
7. Patients with abnormal blood pressure.
8. Patients who cannot freely express their wishes.
9. Patients with sensory impairment due to high degree of peripheral circulatory disturbance caused by diabetes.

Precautions

It is forbidden to use in the following situations:

1. Patients with malignant tumors;
2. Patients with acute suppurative inflammation;
3. Near the heart of the human body;
4. Patients with bleeding tendency;
5. Pacemaker users;
6. Patients with metal implants;
7. Pregnant women and parturients;

>Do not use near the heart, above the neck, head, mouth, genitals or for skin disorders and do not treat feet at the same time soles of feet. Also, do not press internal organs when using. Failure to do so may result in an accident or physical discomfort.

>Do not use it with other medical devices or more than 2 sets of this product, otherwise it may cause nausea or physical discomfort.

>Do not use it for purposes other than the massage described in the instruction manual, otherwise accidents and malfunctions may occur.

>Do not disassemble, repair or modify this product without authorization, otherwise it may cause malfunction or electric shock.

Precautions

! Attention

1. If you feel unwell due to abnormality of this product, please stop using it immediately.
2. If you feel strange on your body or skin, please consult a doctor and follow the doctor's advice.
3. When there is an abnormal situation that cannot operate, please cut off the power immediately and remove it, otherwise high temperature or short circuit may occur have an accident. At the same time, contact customer service for after-sales consultation.
4. Do not plug or unplug the charging cable with wet hands, otherwise it may cause electric shock or accident.
5. Do not use it for children or people who have no ability to express consciousness, otherwise it will cause accidents or cause physical discomfort.
6. Do not use it in a place with high humidity such as a bathroom, or use the product while taking a bath otherwise you may receive a strong electric shock and damage it product.

8. Do not use the product while driving, otherwise, physical discomfort may occur and thereby cause traffic accidents.
9. Do not apply the product to the same site for too long (over 30 minutes), otherwise, muscle fatigue may occur and thereby cause physical discomfort.
10. Please wear the product correctly and symmetrically, and apply it to the cervical sweater, otherwise accidents or physical discomfort may occur.
11. Do not pull the charging cable to unplugging it, otherwise, an electric shock may occur and the product may be damaged.
12. The surface of the product will heat up after operating for a long time, which is a normal phenomenon.

Precautions

! Attention

- Please confirm the button on the main frameworks before using the massager.
- The battery power of the product may be insufficient when being used for the first time or after being left idle for a long time. In this case, please recharge the product according to the instruction to ensure normal operation.
- Please dispose of the product according to the regulations on environmental protection of the city, otherwise, it may cause environmental pollution.
- The product is recommended for people who are 16 years old or above,

Do not use or store the product in the following places:

- 1 Places with direct sunlight
- 2 Places with high temperature and humidity
- 3 Places with water Dusty places
- 4 Places near fireworks
- 5 Places vulnerable to vibration and shock
- 6 Places with strong electromagnetic interference

Tips

- The product will power off automatically after running for 15 minutes. You are recommended not to use the product for over 30 minutes a day and use it as required every week.
- Please place the product at ventilated and dry places after use.
- Do not put other heavy objects on the product when using it.
- Do not pull the rubber plug at the charging port hard, otherwise, it would be easily broken and fall off.

Q&A

If you encounter any problem when using the product, please refer to the following tips:

- > Turn off the power if the product remains in the faulty state or could not work.
- > Users should not disassemble or repair the product by themselves. Any damages caused by such behavior would invalidate the product warranty, and we will not bear any subsequent consequences.

Failure	Possible cause / countermeasure
The massager could not be started	<ul style="list-style-type: none"> • Please check whether the battery power is low. • Please check whether the button has been pressed (the green light would stay on if the massager is powered on).
The massager could not be recharged or the charging light does not work	Please check whether the cable is connected to the mainframe and the power adapter of the mobile phone.
The massager stopped operation suddenly	<ul style="list-style-type: none"> • Please check whether the battery power is insufficient. • The automatic timer (15 minutes) has expired. • If the problem still exists, please stop using the product and contact our maintenance service center.

Cleaning and Maintenance

Before cleaning the neck massager, please be sure to disconnect the power supply, and do not touch the plug with wet hands.

1. Prohibit close to fire, oil, corrosive liquid and various adverse environments.
2. Only the skin-adhesive outer cloth cover can be disassembled and cleaned separately, and parts such as massage components and control panels are prohibited from getting wet.
3. Before cleaning, the outer cloth cover must be disassembled as required to ensure that the outer cloth cover is completely separated from the circuit and other non-cleanable parts.
4. The outer cloth cover is not machine washable, but can be gently hand washed, and the water temperature should not exceed 40 degrees. It is not suitable for drying and drying, but hang to dry naturally.
5. After the outer cloth cover is cleaned, other parts must be installed after the outer cloth cover is completely dry, and then turned on for use.
6. Please keep the product out of the reach of children.
7. Do not store the product in a place exposed to direct sunlight, high temperature, or humidity.
8. Please store the product in a dry and ventilated place.
9. Do not disassemble, repair or modify the product without permission, as it may cause malfunction.

Product Disposal Precautions

Please offer your support for the recycling of rechargeable batteries (lithium-ion batteries)

⚠ Dangers

1. The battery used in this product is a special rechargeable battery and cannot be used in other products. Do not charge the rechargeable battery that has been taken out.

2. Do not throw it into the fire to make it overheated.
3. Do not pierce with nails, impact by gravity or disassemble and modify.
4. Do not connect the positive and negative poles of the battery directly with metal.
5. Do not store it together with necklaces, hairpins, etc.
6. Do not charge, use or place the device in high temperature places such as near the stove or under the hot sun.
7. The above situations can cause heat, fire, rupture and other phenomena.

⚠ Warnings

If you encounter any problem when using the product, please refer to the following tips:

1. Please put the removed rechargeable battery out of the reach of infants and young children to prevent it from being swallowed by mistake.
2. If swallowed by mistake, please seek medical treatment immediately.
3. When the battery liquid leaks, do not touch the liquid with bare hands, and handle it as follows:
 - (1) When the battery fluid enters the eyes, there is a danger of blindness. Do not rub the eyes directly with your hands, and immediately wash them with clean water. Wash and seek medical treatment;
 - (2) When the battery fluid gets on the body or clothes, it may cause skin inflammation or injury, please rinse with clean water, if there is still unsatisfactory. If appropriate, seek medical attention immediately.
 - (3) Only the product can be decomposed when it is discarded. Failure to comply may cause fire, electric shock, injury and other accidents. If repairs are required (except for battery replacement), please consult the store where you originally purchased the product.

.....
The ion battery used in this product is a valuable resource that can be recycled. Please dispose of the ion battery after use. Take it out and hand it over to a professional organization that can recycle and recycle rechargeable batteries.

List of names and contents of toxic and hazardous substances in products

Hazardous Substances	Name of Parts				
	Hardware	PCBboard	Plasticparts	Sponge	Cloth
Lead PbPCB	<input type="radio"/>	<input checked="" type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>
Mercury Hg	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>
Cadmium Cd	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>
Cd Hexavalent chromium Cr(VI)	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>
Polybrominated biphenyls PBB	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>
Polybrominated biphenyls ether PBDE	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>

- : Indicates that the content of this hazardous substance in all homogeneous materials of this part is below the limit requirement specified in GB/T26572.
- : Indicates that the content of the hazardous substance in at least one homogeneous material of the part exceeds the limit requirement specified in GB/T26572.

Environmentally friendly use period of 10 years:

The "environmental protection use period" of this product is 10 years, and the company formulates this period according to the product design characteristics.

Only when using this product under the normal conditions described in the instruction manual of this product "Environmentally friendly use period" is only valid.



Notes

- (1) Parts and their homogeneous materials not listed but contained in the product do not contain any hazardous substances.
- (2) All the parts and their homogeneous materials containing hazardous substances listed in the above table all conform to the strict requirements for the limit on usage of hazardous substances as specified in the Restriction of Hazardous Substances Directive (RoHS) of China. Your safety is assured while using the product.

